


# ISOSCAN2

FDXA/FDXB/HDX

Complies with ISO 11784/85 & 14223

## BASIC user manual

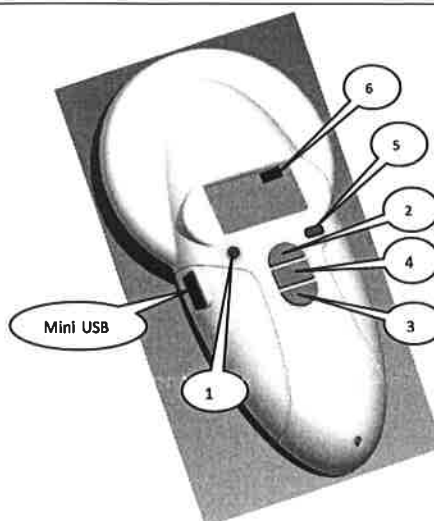

Congratulations, you have just acquired your ISOSCAN2 reader. This reader can read all FDX-B type electronic chips (complying with standard ISO 11784), FDX A chips and HDX chips.

As you will be able to see for yourself, it is extremely simple to use.

### Description of the reader

The reader has four buttons.

- **One round button to switch the device on or off,**

- **Three central buttons.**

The upper button (2) and lower button (3) scroll up and scroll down) are used to scroll the various functions of the reader on the display.

The central rectangular "Validate" button (4) is used to validate the options selected by the user when navigating through the various menus.

At the top left of the device you will see a mini-USB type connector used for:

- + recharging the device via a USB cable,
- + transferring any data stored by the reader.

- **An OLED display comprising four lines of 16 characters.**

### Configuring the reader

After switching on the reader by pressing the On/Off button (1), you will see the message "Read" with an arrow above it.

Press button 2 to display the "Language" menu. The reader is initially configured to operate in English.

To change the language, press button 3 and scroll down the languages until the one you wish to use (for example English) is displayed. You must then validate your choice (button 4).

The following message is displayed for 4 seconds:

English  
Ok

and then the reader displays:

SCAN

The languages available are: English, French, Spanish, Portuguese, Italian, German and Polish.

Once you have selected the language you can switch off the reader by holding down the "On/Off" button.

### Activating the reader

When the reader is switched off, you can switch it on by pressing button 1. The display shows:

SCAN

Press the "Validate" button (4) to activate chip search for 20 seconds.

### Reading

The operator must hold the reader close to the place where the chip is assumed to be, scanning over the area slowly.

Two cases may arise:

#### First case:

A transponder is detected. The reader emits a "beep" and displays the identification number in the following form:

FDXB  
932 274877906744

In this case the transponder is an ISOFDXB type transponder. The identification number takes the form of a series of 15 digits (FDX B) or 10 digits (FDXA). The number remains displayed for two minutes before the reader switches itself off.

**Note: switching off the reader causes the number displayed to be erased from the memory.**

#### Second case:

No transponder is detected. If no transponder has been detected after 20 seconds, the reader emits three "beeps" and displays:

No transponder  
found !

This message will be erased after approximately 15 seconds. The reader will then display:

SCAN

It will switch itself off automatically after two minutes to save battery power.

**In this case, it is recommended to make two or three attempts to read a chip, scanning more widely over the area in which it is assumed to be implanted. Each successive pressing of the button (4) will restart the reading in accordance with the process described above.**

### Recharging the battery

The reader is powered by a Lithium/Ion battery. This battery is designed to allow several thousand reading operations. It can be recharged by connecting the reader to a USB port. A fixed red LED (5) indicates that recharging is in progress and green (5) when it is fully loaded.

The service life of the battery depends not only on its use but also on the environment in which the reader is stored.

An Indicator on the right of the display (6) indicates the battery charge level. When the charge is too low the reader displays the following message:

### Low battery

Only around ten more reading operations will be possible before the reader switches itself off permanently.

To replace the battery, open the battery compartment on the underside of the reader. It is essential to use original batteries only. (Consult your distributor).

### Specifications

Complying with standards ISO 11784/85 and 14223

Dimensions: L 15 cm, W 8 cm, Thickness 3.5 cm

Weight: 150 g

Power supply: 3.7V 1400mA rechargeable battery

2 x 12mm FDXA glass tag scanning distance: 9cm

2 x 12mm FDXB glass tag scanning distance: 9cm

Reads HDX and EM4102 chips

Supplied with a USB cable

Storage temperature: -10° + 60°

Operating temperature: -5° to +40°

### Options

Bluetooth

### Certifications

CE and FCC

### Guarantee

One year parts and labour with return to factory.