

TMS30 Vet

Veterinary Telemetry System

Focused Care with All Heart



Veterinary Telemetry System

The TMS30 Vet telemetry monitoring system provides real-time monitoring of ECG and SpO₂ for ICU inpatients. Especially suitable for those with heart disease and associated risks to help to effectively reduce medical risks and increase the success rate of treatment.

Together, with the CMS, the whole system helps to effectively minimize the manpower costs related to patient coordination and care, facilitating the establishment of a specialized, information-driven ICU within animal hospitals and boosting the hospital's industry competitiveness.

Easy to wear, Simple to use

With its compact and lightweight design, the TMS30 Vet transmitter can be securely fastened with the professional vet telemetry pouch, preventing damage from patient biting and pulling and other actions.

Professional Vet Telemetry Pouch

- Three sizes to meet different animal body shapes
- Concealed cable storage design, lead wires won't be snagged or chewed
- Electrode fixation points ensure fast and precise placement of ECG leads

Compact & Durable Transmitter



-  **1.54"** color display
Real-time display
-  **IPX7** Waterproof
Easy to clean & disinfect
-  **180h** Battery Life
Supports long-term continuous monitoring
-  **170g** Weight
Lightweight and comfortable
-  **1.5m** Drop resistance
No damage from accidental drops

Unique Shave-free Electrode

- Combined electrode design, no shaving required for ECG monitoring
- Metallic pad with gel-storage electrode increases electrode contact area for stable ECG signals
- Disposable gel-storage electrode for max. 48h continuous monitoring

Accurate in Monitoring, Specialized for ICU

Precise Parameter Monitoring

! ECG or ECG+SpO₂

- Anti low perfusion SpO₂ algorithm
- Multi-lead ECG analysis
- Built-in sensor to detect motion condition

🏠 Comprehensive ECG monitoring

- Monitoring of ST segment and QT/QTc
- Up to 22 types of arrhythmias


ECG SpO₂ Intelligent CrozFusion

- Reliable measurements
- Reduce false alarms

Focused Care for All

- Alarm free on transmitter, brings no disturbance to pet
- Audio & visual alarm on CMS to ensure no abnormalities are missed
- Centralized monitoring of multiple patients brings focused care to all patient across departments
- The central station system enables full-disclosure waveform storage, review, and report printing


Multiple Patients' Monitoring

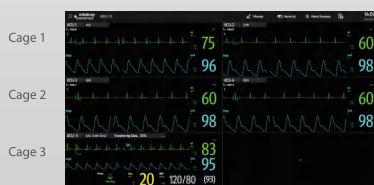

Viewing Telemetry Patient



Full-disclosure Wave Review & Storage

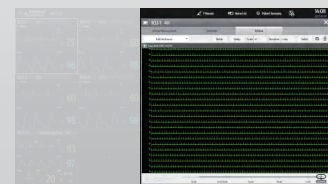


Enhance Professionalism, Help to Increase Revenue

- Seamless care for inpatients with heart disease, brings confidence and satisfaction to pet owner
- Objective and efficient data transmission, reduce human errors
- Enhancing the efficiency of medical staff, provides more focused and comprehensive animal care services

Professional Veterinary ICU Solution

Intensive Care for **Cardiac Hospitalization Patient**

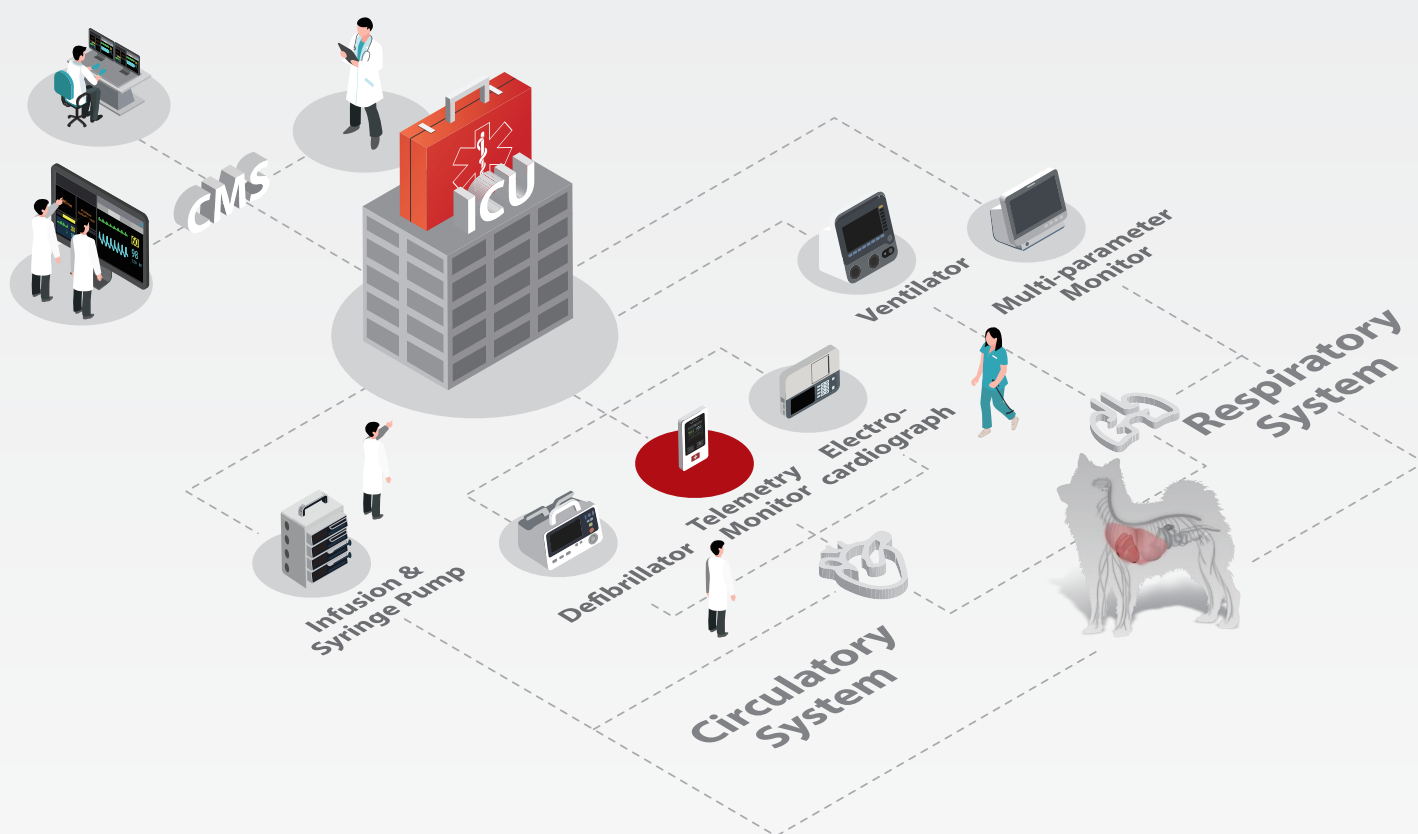

LinkedIn



Facebook



Instagram



Search Mindray Animal Medical

www.mindrayanimal.com
service@mindrayanimal.com

P/N:ENG-TMS30 Vet-210285X4P-20240512
©2022 Shenzhen Mindray Animal Medical Technology Co., LTD. All rights reserved.

 **mindray**
animal care